

## USER GUIDE

### Luxury paraffin bath, ref. FINEX (WK-E004)

Please read this document carefully before using the device. Keep it in an accessible place so that you can read it at any time.

#### General information


The paraffin bath is popular for beauty and home treatments because of its outstanding effect on dry and aged skins, but also on anti-wrinkle treatments.

This easy-to-use equipment is popular among the beauty centres all over the world for the following features: a temperature controller and display, a quick melting, a timing controller and display, a temperature switch in Fahrenheit and Celsius degrees and a large capacity.



1. Lid

2. General display

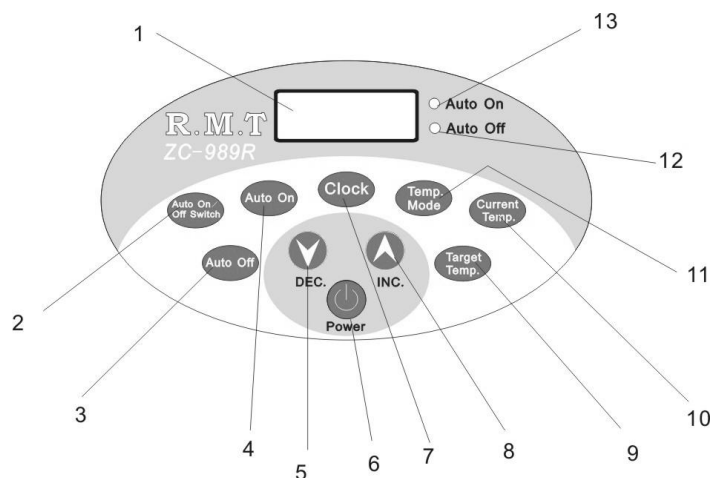

1. Display

2. *Auto on/off switch* : Automatic mode switch
3. *Auto off*: Time setting to end the automatic mode
4. *Auto on*: Time setting to start the automatic mode
5. *DEC.*: Temperature/time down
6. Power button
7. Clock
8. *INC.*: Temperature/time up
9. Target temperature
10. Current temperature
11. Temperature unit setting
12. Automatic mode off indicator
13. Automatic mode on indicator

## Operation

1. To switch the equipment on, press the button “Power”.
2. The first time you use the equipment, the clock needs to be adjusted on time. To do so, press the button “7” and the first two numbers of the clock will shine. Then, by pressing the buttons “5” and “8”, you will be able to adjust the hour. If you press the button “7” again, the last two numbers of the clock will shine. Then, press the buttons “5” and “8” to adjust the minutes. To complete the clock setting, press the button “7” again.
3. To adjust the working mode, press the “2” button. Then, the “Auto on” and “Auto off” indicator will shine. If you press the “4” button, the numbers of the display will shine and you will be able to adjust the time at which the automatic mode has to be switched on by pressing the “5” and “8” buttons. If you press the “3” button, you will be able to set the time at which the automatic mode has to be switched off. If the indicators (“12” and “13”) do not shine, it means that the machine is working in the manual mode.
4. Temperature unit setting. Press the “11” button to select the temperature unit. Then, the “5” and “8” will allow you to choose either the Fahrenheit (F) or the Celsius (C) scale.
5. Target temperature setting. Press the “9” button to set the target temperature. The “5” and “8” buttons will allow you to reach the temperature needed.
6. Current temperature display. Press the “10” button to see the current temperature.

### **Remarks and recommendations**

1. After turning on the automatic mode, do not unplug the equipment. Make sure that the clock is on time.
2. If the automatic mode is not used, make sure that you switch the equipment off after use.
3. The machine can save the temperature that was set last time.
4. Before starting the treatment, check the temperature on your hand to see if it is the proper one and make sure it will not be uncomfortable for the customer.
5. If the skin is damaged, do not apply this treatment.
6. After finishing the treatment, do not pour the remaining paraffin when it is still liquid. Turn off the unit and wait for the paraffin to get colder and become completely solid. Then, turn the equipment on again and wait for the paraffin at the bottom to be melted and you will be able to remove the paraffin block easily.
7. In order to shorten the melting time, choose the maximum temperature. Once the paraffin is completely melted, then adjust the temperature at the proper one.
8. If the power cord is damaged, it should be replaced by one provided by the manufacturer or the authorised technical service.

### **Technical specifications**

Complete melting time: 90 minutes

Adjustable temperature range: 46-60°C

Capacity: 5000ml

Rated power: 200W

Rated voltage: AC 220V ~ 240V